

Cumulative Release Notes for DSW.

Version 4.00

This update covers the changes required for 1st July 2008.

The Diagnosis page now includes a Condition Onset flag. Note that these radio buttons do not represent the same thing in every state. Some states define the flag as “Condition Present on Admission”, others as “Condition with onset during the episode of admitted patient care”. The former means when the flag is set the disease did not occur while the patient was admitted, the later means it did occur while the patient was admitted.

The 3M grouper Interface has been updated to support version 2.4.x of the software, this is the ONLY version of the Grouper currently supported

The test for completeness invoked by the Completed button has been enhanced to check for the Clinician on Admission and Discharge’s MBRN and the Clinician on Admission’s Specialty.

This release supports Overnight Stay. If the hospital has Overnight Stay patients this feature is enabled by the “Overnight” Property in Settings having a Value of “Y”

Installation.

1. Rename DSW.mdb, DSW Data.mde and DSW Episodes.mde to Old DSW.mdb, Old DSW Data.mde and Old DSW Episodes.mde respectively.
2. Copy the new DSW.mdb, DSW Data.mde DSW Episode.mde and DSW Updater.mdb into the same folder.
3. Run the Updater to copy the data from the Old DSW Episodes.mde database into the new version.
4. Start DSW and Relink.

Version 4.01

Features.

1. Support for 3M Core Grouping Software Version 1.1

Installation procedure.

Replace DSW.mdb and relink. Refer to the Separate document “3M Core Grouper 1_1 DS Notes.doc” for instructions on installing / updating the Grouper software and configure it and Day Surgery.

Modified File List

File Name	Version	Modified from last release
DSW.mde	4.01	Yes
DSW Data.mde	2.00	No
DSW Episodes.mde	3.00	No
DSW Updater.mde		No / Not needed

Test Instructions

Test with data files in both “product” and other (e.g. d:\HospData\) folders.

Area of product affected.

3M Grouper Interface

Reference

3M Core Grouper 1_1 DS Notes.doc

Version 4.02

Features.

Adds enhanced support for Overnight Stay patients as follows:

1. An additional field "Same Day Status" has been added to the Discharge page of the Episodes screen. For new episodes this field defaults to Same Day (1).
2. The functionality of the HC21 form has been enhanced. The box on the right hand side of the page labelled Same Day Patients Only will not have any items filled in if the same day status is not "1" i.e. Same Day.
3. The HCP based extracts have been enhanced to report the actual same day status rather than a default value of "1".
4. The HMDS extracts have been enhanced to report a value for Intended Length of Stay based on the value entered for Same Day Status.

Installation procedure.

As changes have been made to DSW Episodes.mde it will be necessary to update existing site's data into the new database. Once this is done replace DSW.mdb and DSW Data.mde, relink if required. Check that all site-specific settings have been preserved.

If the site is doing overnight stay check that their settings table has a Parameter "Overnight" and its value is set to "Y". The effect of having this parameter set to "N" or not having the parameter is to cause the Separation Date to track the Admission Date. A value for the Separation Date different to the Admission Date may still be entered

Modified File List

File Name	Version	Modified from last release
DSW.mde	4.02	Yes
DSW Data.mde	2.01	Yes
DSW Episodes.mde	3.01	Yes
DSW Updater.mde		No

Test Instructions

1. Check that the correct Same Day Status is sent in the HCP extracts.
2. Check that the correct derived value is sent as Intended Length of Stay in the HMDS Extract.
3. Check that the HC21 Claim form is working as per spec.
4. Check that the Updater works correctly

Area of product affected.

1. HC21 Claim Form
2. HCP based data extracts (PHDB, Hospital to Fund)
3. HMDS Extract

Reference

HCP data extract file format

HMDS data extract specification and manual

Version 4.03

Features.

This release corrects a problem with HMDS extracts.

The Clinician at Admission was being used as the Clinician Performing Procedure for all procedures reported including Anaesthetic procedures which is very likely to be incorrect.

As the reporting of Clinician Performing Procedure is optional this release corrects the issue by not reporting the Clinician Performing Procedure for all secondary procedures.

Installation procedure.

From Version 4.02

Replace DSW.mde and Relink.

From versions older than 4.02

Carry out update as per instruction for 4.02 **using the Updater supplied with this release.**

Modified File List

File Name	Version	Modified from last release
DSW.mde	4.03	Yes
DSW Data.mde	2.01	No
DSW Episodes.mde	3.01	No
DSW Updater.mde	2.03	Yes

Test Instructions

1. Check that no Clinician MBRN is reported in all but the Principal Procedure
2. Check that the Updater correctly copies ALL data items. The previous version failed to copy the table of secondary procedures.

Area of product affected.

1. HMDS extracts for 2007 and 2008
2. Updating from previous versions

Reference

HCP data extract file format

HMDS data extract specification and manual

Version 4.04

Features.

Updates data tables for HMDS for 2009

Updates the Multiselect function so Cancel aborts to same functionality as if not enabled.

Updates tge 3M Grouper Interface to support 3M Core Grouping Software V1.3 and AR-DRG V6.0.

Installation procedure.

From Version 4.02, 4.03

Replace DSW.mdb and DSW Data.mde and Relink.

From versions older than 4.02

Carry out update as per instruction for 4.02.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	4.04	Yes
DSW Data.mde	2.02	Yes
DSW Episodes.mde	3.01	No
DSW Updater.mde	2.03	No / Not required

Area of product affected.

1. HMDS COB, Clinician Speciality data tables
2. 3M Core Grouping Software
3. Multiselect.

Reference

3M Core Grouping Software V1.4 manual

HMDS data extract specification and manual.

HMDS Update documentation for 01/07/09

Version 5.00

Features.

Updates data tables for ICD-10-AM Version 7

Update HCP based extracts for HCP version 8.00.

Updates the 3M Grouper Interface to support 3M Grouping Software Version 2010.2 and ICD-10-AM Version 7

Installation procedure.

From Version 4.02, 4.03

Replace DSW.mdb and DSW Data.mde and Relink.

Either update the Local Health Funds and Group Funds manually or by using the generic fund updater. If the update is undertaken after 1st July 2010 then the generic updater is recommended as it will update completed episodes to the new fund codes.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	4.04	Yes
DSW Data.mde	2.02	Yes
DSW Episodes.mde	3.01	No
DSW Updater.mde	N.A.	No. Use generic updatere

Area of product affected.

1. ICD-10-AM coding tables
2. 3M Grouping Software interface
3. HCP based extracts

Reference

Appendix_09_HMDS_Interface_File 2010.doc

HMDS_July_2010_Changes_Private Hospitals_v1.0.doc

HCP data specifications _hospital to insurer_ 2010-11.pdf

Version 6.00

Features.

- Updates HC21 forms and HCP extracts to handle “ADA” item codes and show correct age.
- Updates the 3M Grouper Interface to support 3M Grouping Software Version 2010.3
- Collects and reports Baby Weight and Unplanned Return to Theatre for HMDS
- Updates MBRN for 13 character nation registraion numbers

Installation procedure.

From Version 5.00

As changes have been made to DSW Episodes.mde it will be necessary to update existing site's data into the new database using the supplied Updater. Once this is done replace DSW.mdb and DSW Data.mde, relink if required. Check that all site-specific settings have been preserved

As this fund is also part of AHSA the new code must also be added to the Grouped Health Funds. Do this from Advanced>Edit Datasets>Grouped Health Funds.

Work to be done by the hospital

As HMDS have changed from using state based Medical Board Registration Numbers to the national registration number the dataset needs to be updated. This is done through Advanced>Edit Datasets>Attending Doctors. Navigate through each doctor and replace the MRBN with their national registration number. The form in question is shown below.

Attending Doctor Details

Name:

Specialty:

Address 1:

Address 2:

Suburb: Post Code: State:

Phone Number:

Provider Number:

MBRN:

Valid: ☒

OK

Record: 1 of 1

Use these navigation buttons to browse through every doctor.

WARNING: Changes to MBRN numbers must not be done until after all episodes prior to 01/07/11 have been exported to HMDS.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	6.00	Yes
DSW Data.mde	4.01	Yes
DSW Episodes.mde	2.02	Yes
DSW Updater.mde	3.00	Yes

Area of product affected.

1. 3M Grouping Software interface
2. HCP based extracts.
3. HMDS Extracts

References

HMDS Manual 2010.

HCP extract specification 0700

Version 6.01

Features.


This release addresses the following issue:

The PHIC in circular PHI 62/11 have mandated that episodes for the following health funds are reported with fund code "BPU" from 1st October 2011 :

- Hospital Benefits Association Ltd (HBA)
- Mutual Community Ltd (MCL)
- MBF Australia (MBF)
- MBF Alliances Pty Ltd (SGI)

Installation procedure.

1. Replace the existing DSW Data.mde file with the newest version.
This version has a new fund in it with code "BUP" and description "Bupa Australia".
This can be done at any time prior to 1st October 2011.
2. Edit the Local Health Fund table as follows after close of business on 30th September 2011 and before opening for business :
 - a. Open the Local Health Fund table through Advanced > Edit Datasets > Local Health Funds.
 - b. You should now see a screen similar to the following:

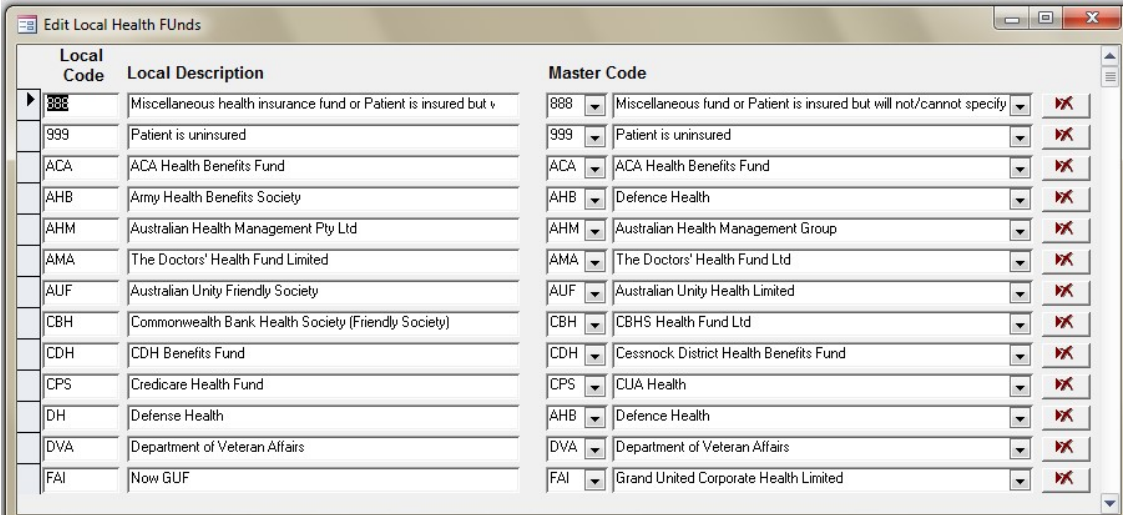

Local Code	Local Description	Master Code	Master Description
888	Miscellaneous health insurance fund or Patient is insured but v	888	Miscellaneous fund or Patient is insured but will not/cannot specify
999	Patient is uninsured	999	Patient is uninsured
ACA	ACA Health Benefits Fund	ACA	ACA Health Benefits Fund
AHB	Army Health Benefits Society	AHB	Defence Health
AHM	Australian Health Management Pty Ltd	AHM	Australian Health Management Group
AMA	The Doctors' Health Fund Limited	AMA	The Doctors' Health Fund Ltd
AUF	Australian Unity Friendly Society	AUF	Australian Unity Health Limited
CBH	Commonwealth Bank Health Society (Friendly Society)	CBH	CBHS Health Fund Ltd
CDH	CDH Benefits Fund	CDH	Cessnock District Health Benefits Fund
CPS	Credicare Health Fund	CPS	CUA Health
DH	Defense Health	AHB	Defence Health
DVA	Department of Veteran Affairs	DVA	Department of Veteran Affairs
FAI	Now GUF	FAI	Grand United Corporate Health Limited

c. Scroll down the table to find the entry for Hospital Benefits Association under Local Code and/or Local Description.

d. The entry for HBA is highlighted below:

Edit Local Health Funds			
Local Code	Local Description	Master Code	
CBH	Commonwealth Bank Health Society (Friendly Society)	CBH	CBHS Health Fund Ltd
CDH	CDH Benefits Fund	CDH	Cessnock District Health Benefits Fund
CPS	Credicare Health Fund	CPS	CUA Health
DH	Defense Health	AHB	Defence Health
DVA	Department of Veteran Affairs	DVA	Department of Veteran Affairs
FAI	Now GUF	FAI	Grand United Corporate Health Limited
GEH	Now AHM	GEH	
GMF	Goldfields Medical Fund (Inc)	GMF	Goldfields Medical Fund (Ceased Trading)
GMH	Geelong Medical & Hospital Benefits Association Ltd	GMH	Geelong Medical & Hospital Benefits Association Ltd
GUF	Grand United Friendly Society	GUF	Grand United Friendly Society (Ceased Trading)
HBA	Hospital Benefits Association Ltd	HBA	BUPA Australia Health Pty Ltd - Discontinued 01/10/11
HBH	Hospital Benefits Fund of Western Australia (Inc), The	HBH	HBH Health Funds Inc - Discontinued 01/10/11
HCF	Hospitals Contribution Fund of Australia Ltd, The	HCF	Hospitals Contribution Fund of Australia Ltd, The

e. Select the “BUP” Master Code clicking the dropdown arrow as shown below:

CBH	Commonwealth Bank Health Society (Friendly Society)	CBH	CBHS Health Fund Ltd	
CDH	CDH Benefits Fund	CDH	Cessnock District Health Benefits Fund	
CPS	Credicare Health Fund	CPS	CUA Health	
DH	Defense Health	AHB	Defence Health	
DVA	Department of Veteran Affairs	DVA	Department of Veteran Affairs	
FAI	Now GUF	FAI	Grand United Corporate Health Limited	
GEH	Now AHM	GEH		
GMF	Goldfields Medical Fund (Inc)	GMF	Goldfields Medical Fund (Ceased Trading)	
GMH	Geelong Medical & Hospital Benefits Association Ltd	GMH	Geelong Medical & Hospital Benefits Association Ltd	
GUF	Grand United Friendly Society	GUF	Grand United Friendly Society (Ceased Trading)	
HBA	Hospital Benefits Association Ltd	HBA	BUPA Australia Health Pty Ltd - Discontinued 01/10/11	
HBH	Hospital Benefits Fund of Western Australia (Inc), The	HBH	HBH Health Funds Inc - Discontinued 01/10/11	
HCF	Hospitals Contribution Fund of Australia Ltd, The	HCF	Hospitals Contribution Fund of Australia Ltd, The	

f. The entry should now appear as follows:

HBA	Hospital Benefits Association Ltd	BUP	Bupa Australia
-----	-----------------------------------	-----	----------------

g. Repeat for the remaining 3 health funds listed above.

- After completing the above changes check that when each of the above health funds are selected in the Claims tab of the episode record that the Master Fund Code is displayed as BUP as shown below:

The screenshot shows a software window with several tabs: Demographics, Admission, Discharge, Claim, Fees, Diagnoses, Ext Causes, and Proceed. The 'Claim' tab is active. In the 'Local Fund' section, there is a dropdown menu showing 'HBA' and 'Hospital Benefits Association Ltd'. To the right of this dropdown is a button labeled 'BUP'. A blue box labeled 'Master Fund Code' has an arrow pointing to the 'BUP' button. Below the 'Local Fund' section are several other fields: Product code, Membership #, Fund cover, Member relationship, Member first name, Member surname, Member address, Member suburb, Member home phone, and Member work phone.

Note that only new episodes will be affected. It is for this reason that the changes must be scheduled after 30th September 2011 and before opening on or after 1st October 2011.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	6.01	No
DSW Data.mde	4.01	Yes
DSW Episodes.mde	2.02	Yes
DSW Updater.mdb	3.00	N/A

Area of product affected.

- HCP based extracts.

References

HIC Circular PHI 62/11

Version 6.02

Features.

Corrects a problem with the HMDS extracts post 30/06/11 which had I data item being placed at the wrong offset in the "SUM" inliner record.

Installation procedure.

From Version 6.01

Replace DSW.mdb and relink.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	6.02	Yes
DSW Data.mde	4.02	No
DSW Episodes.mde	2.03	No
DSW Updater.mde	3.00	No / N.A.

Area of product affected.

1. HMDS Data Extracts

References

HMDS Manual December 2010.

Version 7.00

Features.

1. 3M Grouper interface updated for 3M CGS Version 2012.2
2. HCP 09.00 implemented for extracts after 31/10/12

Installation procedure.

From Version 6.01 / 2

Replace DSW.mdb, DSW Data.mde and relink.

Update Local Health Funds and Grouped Health Funds as needed to accommodate new health fund HEA, Check that CWH is included in Local Health Funds.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	7.00	Yes
DSW Data.mde	4.03	Yes
DSW Episodes.mde	2.04	Yes, new installs only
DSW Updater.mde	3.00	No / N.A.

Area of product affected.

1. HCP Based extracts
2. 3M Core Grouping Software interface

References

HMDS Manual December 2010

Installing 3M Grouper Version 2012.doc

HCP 09.00 Extract specification.

Version 7.01

Features.

3. Input validation added to Admission Weight
4. Catch All code added for HCP 2013 extracts for Neonate Weight and Mech Ventilation Hours

Installation procedure.

From Version 7.01

Replace DSW.mdb and relink.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	7.01	Yes
DSW Data.mde	4.03	No
DSW Episodes.mde	2.04	No
DSW Updater.mde	3.00	No / N.A.

Area of product affected.

1. HCP Based extracts
2. Episodes data input form

References

HCP 09.00 Extract specification.

Version 8.00

Features.

5. ICD-10-AM tables updated to Version 8
6. 3M CGS Interface updated for Version 2012.3 and DRG Version 7.0
7. Abortion Notification System Extract capability added
8. HCP based extracts updated for new version of ICD-10-AM

Installation procedure.

- Make backups
- Rename DSW Episodes.mde to Old DSW Episodes.mde
- Copy the DSW Updater into the same folder and run it to import records from the Old DSW Episodes file into the new.
- Replace DSW.mdb and DSW Episodes.mde with the new versions
- Start DWE.mdb by double clicking it. Ignore any error messagea and Relink the two backend databases using Advanced>Relink Databases

Modified File List

File Name	Version	Modified from last release
DSW.mdb	8.00	Yes
DSW Data.mde	5.00	Yes
DSW Episodes.mde	3.00	Yes
DSW Updater.mde	8.00	Yes

Version 8.01

Features.

9. HMDS “SUM” inliner record changed to be the correct length

Installation procedure.

- Replace and Relink DSW.mdb

Modified File List

File Name	Version	Modified from last release
DSW.mdb	8.01	Yes
DSW Data.mde	5.00	No
DSW Episodes.mde	3.00	No
DSW Updater.mde	8.00	No

Version 9.00

Features.

10. Medicare IRN is included in HMDS extracts post 30th June 2014
11. All Extracts display a progress meter.
12. DSW will run under Office 64 or 32 bit versions.

Installation procedure.

- Replace and Relink DSW.mdb

Modified File List

File Name	Version	Modified from last release
DSW.mdb	9.00	Yes
DSW Data.mde	5.00	No
DSW Episodes.mde	3.00	No
DSW Updater.mde	8.00	No

Version 9.01

Features.

13. Corrects an issue that can arise with HMDS extracts for 2014 when patient does not report Medicare Card Number and patient is covered by a DVA card.
14. The validation rules for completing an Episode remain unchanged i.e. the patient's fund must be DVA if the Medicare and IRN is blank.

Installation procedure.

- Replace and Relink DSW.mdb

Modified File List

File Name	Version	Modified from last release
DSW.mdb	9.00	Yes
DSW Data.mde	5.00	No
DSW Episodes.mde	3.00	No
DSW Updater.mde	8.00	No

Version 10.00

Features.

- 15. Updates ICD-10-AM tables to Version 9.
- 16. Updates 3M CGS Interface to support version 2015.1 and DRG Version 8.0
- 17. Updates HCP based extracts to reflect latest version of ICD-10-AM.

Installation procedure.

- Replace DSW.mdb and DSW Data.mde and Relink

Modified File List

File Name	Version	Modified from last release
DSW.mdb	10.00	Yes
DSW Data.mde	6.00	No
DSW Episodes.mde	3.00	No
DSW Updater.mde	N.A.	N.A.

Version 11.00

Features.

1. Changes PHDB and HCP file naming as required for upload to the new portal.
2. Removes the test/production and resubmit flags from the PHDM extract header.
3. Allows for Medical Record Numbers of 20 character length.

Installation procedure.

1. Make backups of **everything**.
2. Transfer data in the existing DSW Episodes.mde file to the new version.
 - i. Rename the existing DSA Episodes.mde file to Old DSW Episodes.mde
 - ii. Copy the DSW Updater.mdb file to the system.
 - iii. Run the Updater. Note this version runs much faster checking integrity.
3. Copy the new DSW.mdb file to the system replacing the previous versions.
4. Relink and check all data is still present. Give some attention to data in the Settings of DS.

Modified File List



File Name	Version	Modified from last release
DSW.mdb	11.00	Yes
DSW Data.mde	6.00	No
DSW Episodes.mde	4.00	Yes
DSW Updater.mde	9.00	Yes

Version 11.01

Features.

1. New health fund CBHS Corporate Health Pty Ltd with code “CBC” added to fund table in DS Data.
2. Corresponding entry added into Local Funds in Episodes.
3. Fund code CBC added to list of Group Funds under AHS

Installation procedure.

1. Replace DSW Data.mde with version 6.01.
2. Open DS.
3. Under Advanced>Edit Datasets>Local Health Funds add a new entry for CBC.
A screenshot of a software interface showing a table with two columns. The first column contains the code 'CBC' and the second column contains the name 'CBHS Corporate Health'. There are dropdown arrows and a delete icon (a red 'X') on the right side of the table.
4. Under Advanced>Edit Datasets>Grouped Health Funds add CBC as a member of AHS.
A screenshot of a software interface showing a table with three columns. The first column contains the code 'AHS', the second column contains the code 'CBC', and the third column contains the name 'CBHS Corporate Health'. There are dropdown arrows and a delete icon (a red 'X') on the right side of the table.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	11.01	Yes, only applicable for new installs
DSW Data.mde	6.01	Yes
DSW Episodes.mde	4.00	No
DSW Updater.mde	9.00	No, N.A.

Version 11.02

Features.

1. Adds the new Health Funds ESH and NMW to the list of health funds

Installation procedure.

1. Make backups of **everything**.
2. Replace DSW Data.mde with the new version (6.02).
3. Start DSW and from Advanced>Edit Datasets>Local Health Funds add new entries for the two new health funds:

Local Code	Local Description	Master Code
ESH	Emergency Services Health	ESH Emergency Services Health
NMW	Nurses & Midwives Health	NMW Nurses & Midwives Health

4. Add the two new health funds into the Grouped Health Funds from Advanced>Edit Datasets>Grouped Health Funds as members of AHS:

AHS	ESH	Emergency Services Health
AHS	NMW	Nurses & Midwives Health

Modified File List

File Name	Version	Modified from last release
DSW.mdb	11.02	No
DSW Data.mde	6.02	Yes
DSW Episodes.mde	4.01	Yes - New installs only

Version 12.00

Features.

1. ICD10-AM data tables updated to V10
2. HCP based extracts upgraded to reflect change to ICD10-AM
3. Grouper interface upgraded for AR-DRG V9.0 and 3M CGS Version 2017.1
4. New Health Fund “myOwn” added

Installation procedure.

1. **Backup all DSW files.**
2. Replace DSW Data.mde with version 7.00
3. Replace DSW.mdb with version 12.00
4. Start DSW and from Advanced>Relink relink the Episodes and Data databases.
5. From Advanced>Edit Datasets>Local Health Funds add new entries for the new health fund code MYO, name “myOwn”:

Local Code	Local Description	Master Code	Fund code for PRS/2 exports
MYO	myOwn	MYO myOwn	MYO

6. Add the new health fund into the Grouped Health Funds from Advanced>Edit Datasets>Grouped Health Funds as members of AHS:

Service	Fund	Delete
AHS LHM	Peoplecare Health Insurance	
AHS MIM	Queensland Country Health (Obsolete)	
AHS MYO	myOwn	
AHS NHB	Navy Health Ltd	
AHS NIB	NIB Health Funds Ltd	
AHS NMW	Nurses & Midwives Health	

Modified File List

File Name	Version	Modified from last release
DSW.mdb	12.00	Yes
DSW Data.mde	7.00	Yes
DSW Episodes.mde	4.02	Yes, only applicable for new installs

Version 13.00

Features.

1. ICD10-AM data tables updated to V11
2. HCP based extracts upgraded to reflect change to ICD10-AM
3. PHDB extracts now include patient details

Installation procedure.

1. **Backup all DSW files.**
2. Replace DSW Data.mde with version 8.00
3. Replace DSW.mdb with version 13.00
4. Start DSW and from Advanced>Relink relink the Episodes and Data databases.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	13.00	Yes
DSW Data.mde	8.00	Yes
DSW Episodes.mde	4.02	No

Version 13.02

Features.

1. Grouper interface
2. HCP based extracts upgraded to reflect change to ICD10-AM
3. PHDB extracts now include patient details

Installation procedure.

1. **Backup all DSW files.**
2. Replace DSW.mdb with version 13.02
3. Start DSW and from Advanced>Relink relink the Episodes and Data databases.
4. Assuming the site is running 3M CGS and has updated it to at least 2019.1.x update the Grouper version in Settings.

Note:

The Grouper interface has been updated to support CGS 2019.2 for AR DRG v10 however as it was not available for testing at release of DSW 13.02 it is possible there will be issues.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	13.02	Yes
DSW Data.mde	8.00	Yes
DSW Episodes.mde	4.02	No

Version 13.03

Features.

1. Adds Date Of Birth (DOB) to the find dialog
2. Provides a facility to export a header only file for HCP extracts
3. Correctly flags fund Medibank Private with code "MBP" as obsolete

Installation procedure.

1. **Backup all DSW files.**
2. Replace DSW.mdb with version 13.03
3. Replace DSW Data.mde with version 8.01.
4. Start DSW and Relink

Note:

The HMDS file contains no header or trailer record so Header Only file are not possible.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	13.03	Yes
DSW Data.mde	8.01	Yes
DSW Episodes.mde	4.02	No

Version 13.03.02

Features.

1. Adds revised descriptions for diagnosis codes U07.0, U07.1, U07.2 related to COVID-19.

Installation procedure.

From 13.03.01

1. Backup **everything**.
2. Replace DSW Data.mde with new file
3. Start DSS and under Advanced click About. Check that Data tables version shows 8.02

Note:

This update is optional as all COVID-19 codes are already in the data table. The update only updates their description.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	13.03	No
DSW Data.mde	8.02	Yes
DSW Episodes.mde	4.02	No

Version 13.04

Features.

1. Adds revised descriptions for diagnosis codes U07.0, U07.1, U07.2 related to COVID-19.
2. Supports HCP version 11 for extracts post 01-07-2020

Installation procedure.

From 13.03.01

1. Backup **everything**.
2. Replace DSW Data.mde and DSW.mdb with new files
3. Start DSS and under Advanced click About. Check that Data tables version shows 8.02
4. Under Advanced>Relink select Data and Episodes and to relink

Note:

This update is optional as all COVID-19 codes are already in the data table. The update only updates their description.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	13.04	Yes
DSW Data.mde	8.02	Yes
DSW Episodes.mde	4.02	No

Version 13.05

Features.

1. Addresses an issue with the justification of DRGVer in HCP Version 11 extracts.

Installation procedure.

From 13.04.01

1. Backup **everything**.
2. Replace DSW.mdb with new file
3. Under Advanced>Relink select Episodes and relink

Note:

This update is optional as all COVID-19 codes are already in the data table. The update only updates their description.

Modified File List

File Name	Version	Modified from last release
DSW.mdb	13.04	Yes
DSW Data.mde	8.02	Yes
DSW Episodes.mde	4.02	No

Version 14.00

Features.

1. Updates ICD10-AM tables to V12 in DSW Data.mde
2. Updates all HCP based extracts to reflect those changes
3. For new installs the Group Funds table in DSW Episodes has been updated.

Installation procedure.

From 13.04.01

1. Backup **everything**.
2. Replace DSW.mdb and DSW Data.mde with the new files
3. Under Advanced>Relink relink DSW.Episodes.mde and DSW.Data.mde

Modified File List

File Name	Version	Modified from last release
DSW.mdb	14.00	Yes
DSW Data.mde	9.00	Yes
DSW Episodes.mde	4.03	Yes but only needed for new installs